

Packing Specification

Description: Hurrycane

Product REF: HCANE-BK, HCANE-BL, HCANE-PR, HCANE-RD


Raised By: Darren Kent

Date: 3/10/22

INSTRUCTION MANUAL

Colour, quad-language Instruction Manual to be folded and inserted inside the Retail Pack. Manual comprises of two sheets, each 141 x 610mm. One sheet to have English on one side and Spanish printed on the other side. The second sheet to have French on one side and German on the other. Folded size: 141 x 102mm

File Name: IFU – Hurrycanes (Quad-language) 2022-07 Rev3 UKCA

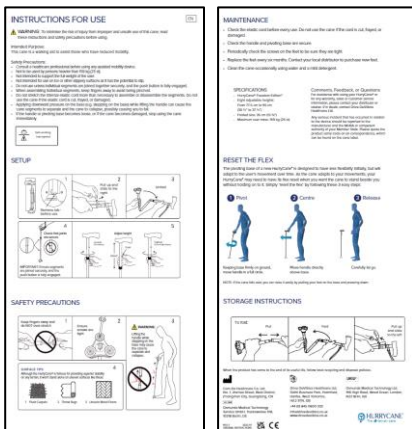

PRODUCT LABEL

Label Size: 80mm x 25mm. Colour: Silver or White.

Serial Number to be added next to SN symbol. Format 18CYMMXXXXXX,

where YY = Year, MM = Month and XXXXXX = unique sequential number, e.g. 18C221000001

Date of Manufacture to be added next to factory symbol. Format YYYY-MM, e.g. 2022-10

Label to be stuck at the bottom of the cane.

File Name: Product Labels – Hurrycane Range (3-10-22)



IMPORTER LABEL

Label Size: 50 x 25mm. Colour: Silver or White.
Label to be on the bottom of the cane, next to the Product Label.


File Name: Importer Label – Drive GB & DE (27-1-22)



RETAIL PACK

Hurricane fixed to colour Retail Pack, with the top of the handle level with the top of the card.
The edges of the pack to be secured with clear sticky pads, in 3 positions; top, side and bottom.

File Name: Retail Pack – Hurricanes (Quad-language) 2022-09 Rev3

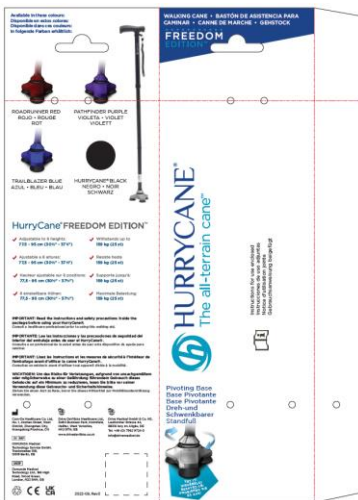

BARCODE LABEL

Label Size: 30 x 30mm, Colour: White,
Label to be stuck to the Retail Pack in position indicated on the artwork.

File Name: Barcode Labels – Hurricane Range (24-1-22)



Barcode label is specific to the colour of the Hurricane

ANTI-RATTLE PROTECTOR

A clear plastic cover to be fitted around the knurled anti-rattle collar, to protect against transit damage.



TUBE PROTECTOR

A clear plastic strip, with 'REMOVE' on it, to be inserted into the end of the folded tube. The strip to be positioned so that it prevents the tube flared end from causing transit damage.

POLYTHENE BAG

The packaged cane to be inserted into a polythene bag, with the bag end folded over and taped closed.

The bag to have a suffocation warning in EN, ES, FR and DE.

WARNING:

TO AVOID DANGER OF SUFFOCATION,
KEEP THIS PLASTIC BAG AWAY
FROM BABIES AND CHILDREN.
DO NOT USE THIS BAG IN CRIBS,
BEDS, CARRIAGES OR PLAYPENS.
THIS IS NOT A TOY.

POLYTHENE BAGS

Note: All polythene bags should preferably contain $\geq 30\%$ recycled content.
A Proof of Specification will be required for the polythene bags used.

CARTON

Plain brown, 5-ply cardboard carton.
Carton Size: 600 x 465 x 390mm (outside dims)

Hurricane canes packed 25 pcs per carton.

Interlocked cardboard dividers to be used to create 5 x 5 spaces.



25 x Serial Number Barcode Labels to be stuck onto the end of the carton



SERIAL NUMBER BARCODE LABELS

Label Size: Approximately 50 x 15mm, Colour: white

The Outer Carton needs to have 25 x Serial Number Barcode Labels stuck on one end of the carton. The Serial Numbers on the 25 labels must match the Serial Numbers shown on the 25 canes that are in each specific Outer Carton.

File Name: Serial Number Barcode Label Template – Com-Da (22-11-21)



Serial Number Format: 18CYMMXXXXXX, where 18C = Factory Code, YY = Year, MM = Month and XXXXX = Unique Sequential Number

SHIP LABELS

Label Size: 140 x 120mm. Colour: White

2x Ship Label to be applied to each carton, on diagonally opposite corners.

Supplier to add the following details to each ship label:

LOT Number, format 18CYYYYMM, e.g. 18C202210

Date of Manufacture, format YYYY-MM, e.g. 2022-10

P.O.#, from Drive purchasing department.

File Name: Ship Labels – Hurrycane Range (3-10-22)

